



Quality from Germany



The **NOVO X 222/333** models represent the next logical step in the expansion of the NOVO X range. They combine exceptional comfort, safety and hygiene as well as all the other first-class features offered by the entire NOVO X range with a patient scale. The elegant display unit can be both swivelled from left to right and fully rotated and is height adjustable – thus guaranteeing optimal user comfort as well as allowing the staff to read the display.

Additional features:

- Weighting range: 0,2 - 200 kg, Verification scale: 0,1 kg
- **New:** with dialyses function
- Flex of the indicator is removable
- Indicator can be run with accumulator



All Likamed products are made in Germany and received certification marks like "TÜV" for top quality. Certification according to ISO 9001 and EN ISO 13485 is the appreciation for innovative and useful developments and for an outstanding product quality.



NOVO X 333

the top model as completely equipped treatment chair which leaves virtually nothing to be desired.

- Three electric drives for a separate adjustment of the back- and the foot part up to the shock and bed position, as well as height adjustment of the complete laying surface.
- Hand switch with „blind symbols“ (only NOVO X 333).

NOVO X 222

with all the features of the type 333 without height adjustment.


- Two electric drives for a separate adjustment of the back- and the foot part.
- Both models are equipped with a flame tested, resistant to disinfectant and toxic tested imitation leather upholstery.

Standard colours:



- upholstery medium grey with matching grey arm and foot rests. Upholstery pastel blue with matching dark blue arm and foot rests. Frame and panelling in grey colour

Other upholstery colours:

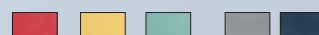

- red, yellow, mint, combined with grey or blue arm and foot rests. Frame and panelling in grey colour (please ask for our special brochure).

Options scale:

- Accumulator, charging cable
- Data interface

For all other options and color samples please see our special brochure.

You will find details of all features in our NOVO X 220/330 brochure.



*



*

* pictures shows NOVO X without scale and without relax-pillow.

Technical specifications and Guarantee

Armrest width	15 cm
Armrest length	57 cm
Laying surface width	60 cm
Laying surface length	213 cm
Total length with extended footrest	220 cm
Sitting - laying height NOVO X 333	56-82 cm
Sitting - laying height NOVO X 222	60 cm
weight	90 kg
Padding and neck support upholstery:	imitation leather, flame retardant
colour: Footpart and armrest:	standard darkblue or darkgrey
Colour upholstery:	lightgrey and pastel blue
frame colour	lightgrey RAL 7047
4 Wheels lockable	Ø = 10 cm
Transformer	protection class II
Transformer rated output primary	230 / 50 Hz
	110 / 50 Hz
	secondary 25 V
fuses	1 AT

on duration	ca. 10%
protection class	II
EMC (electrical safety)	DIN EN 60601/1, 60601/1-2
NOVO X 333	3 x 24 Volt drive with limit circuit breaker
NOVO X 222	2 x 24 Volt drive with limit circuit breaker
Protection category	IP X4 (splash-proof)
Maximum load	175 kg
Ambient temperature	3 - 40°C
Relative air humidity	30 - 75 %
Mains cable rolled	2 m
Mains cable stretched	3,5 m
Tilting angle back section	up to 75°
Tilting angle foot part	up to 35°
Guarantee	2 years
Medical product of class	I
equipotential bonding	
Protection class type CF	optional